


Trilogy Vent

Passive Single Limb Circuit Assembly

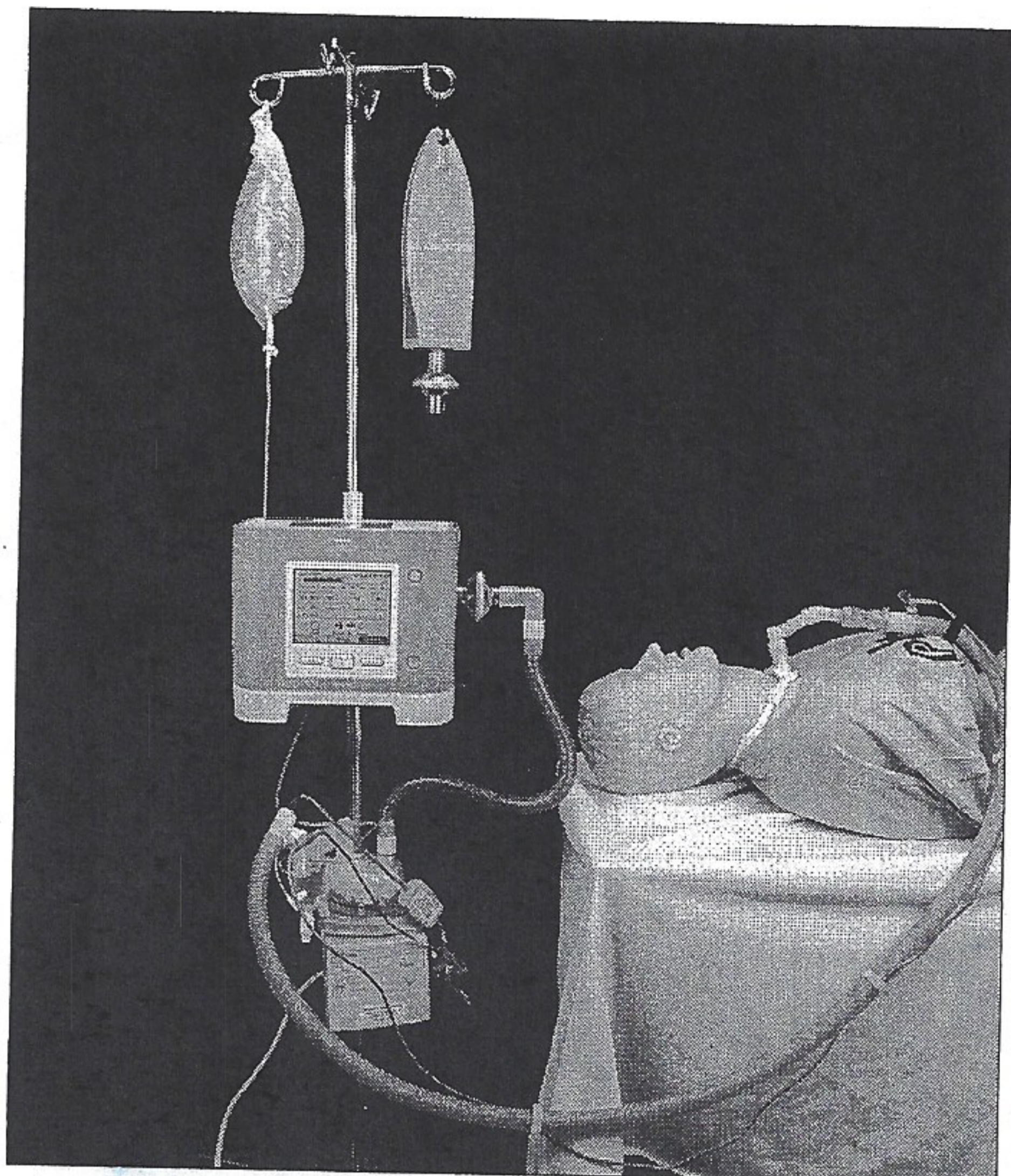

What Does a Trilogy Vent Do?

A Trilogy vent (ventilator) helps your child breathe. You connect the circuit assembly from the vent to your child's trach tube hub. The vent continuously blows air through the circuit assembly and into your child's lungs. A humidifier or heat moisture exchanger (HME) adds warmth and humidity to the air from the vent.


What's Inside	Page
Important Safety Information	2
Equipment	3
Other Equipment Provided Based on Orders and Needs	3-4
Supplies	5-7
Parts of the Trilogy Vent	8
Parts of the Humidifier	8
Attaching Equipment to the IV Pole	9
Setting Up the Humidifier and Water Chamber	10-11
Setting Up the Circuit Assembly	12-14
Connecting the Temp Probe and Pigtail to the Humidifier and Long Tubing Circuit	15
Using the Vent	16-19
Ending Use of the Vent	20
Using the Vent and the HME	21-22
Giving a Nebulizer Treatment While Using the Vent	23-24
Using the Detachable Battery	25
Using the External Battery	26
Power Sources	27
Humidifier Alarms	28
Vent Alarms	29-33
Troubleshooting	34-37
Cleaning—Equipment	38
Cleaning—Supplies	39-40



Read all manufacturer's warnings and cautions before operating equipment for the first time.

Setting Up the Humidifier and Water Chamber


NOTE! Clean your hands before you follow the steps below.

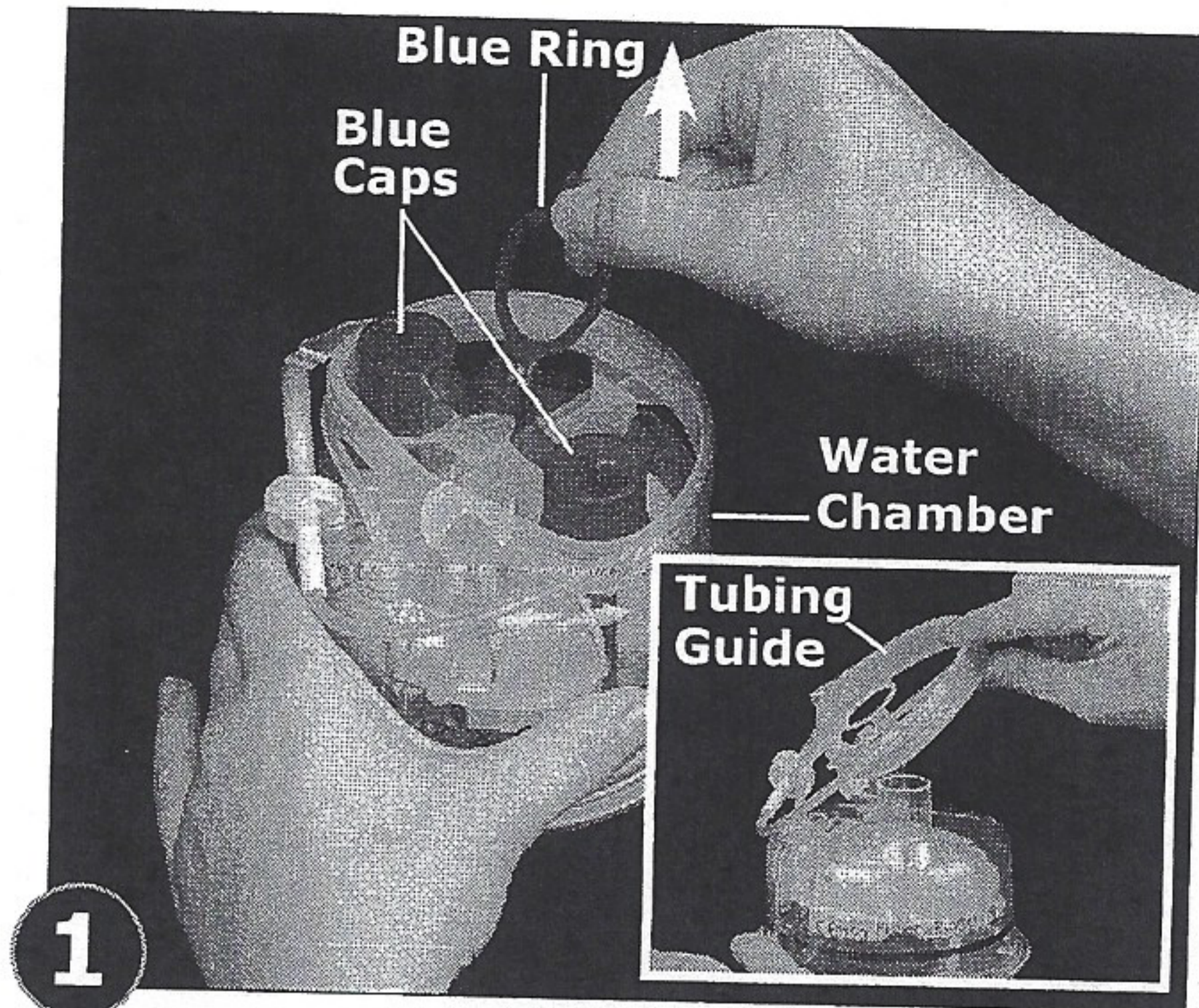

1

Pull up blue ring to remove blue caps from top of water chamber.


Remove tubing guide from top of water chamber.

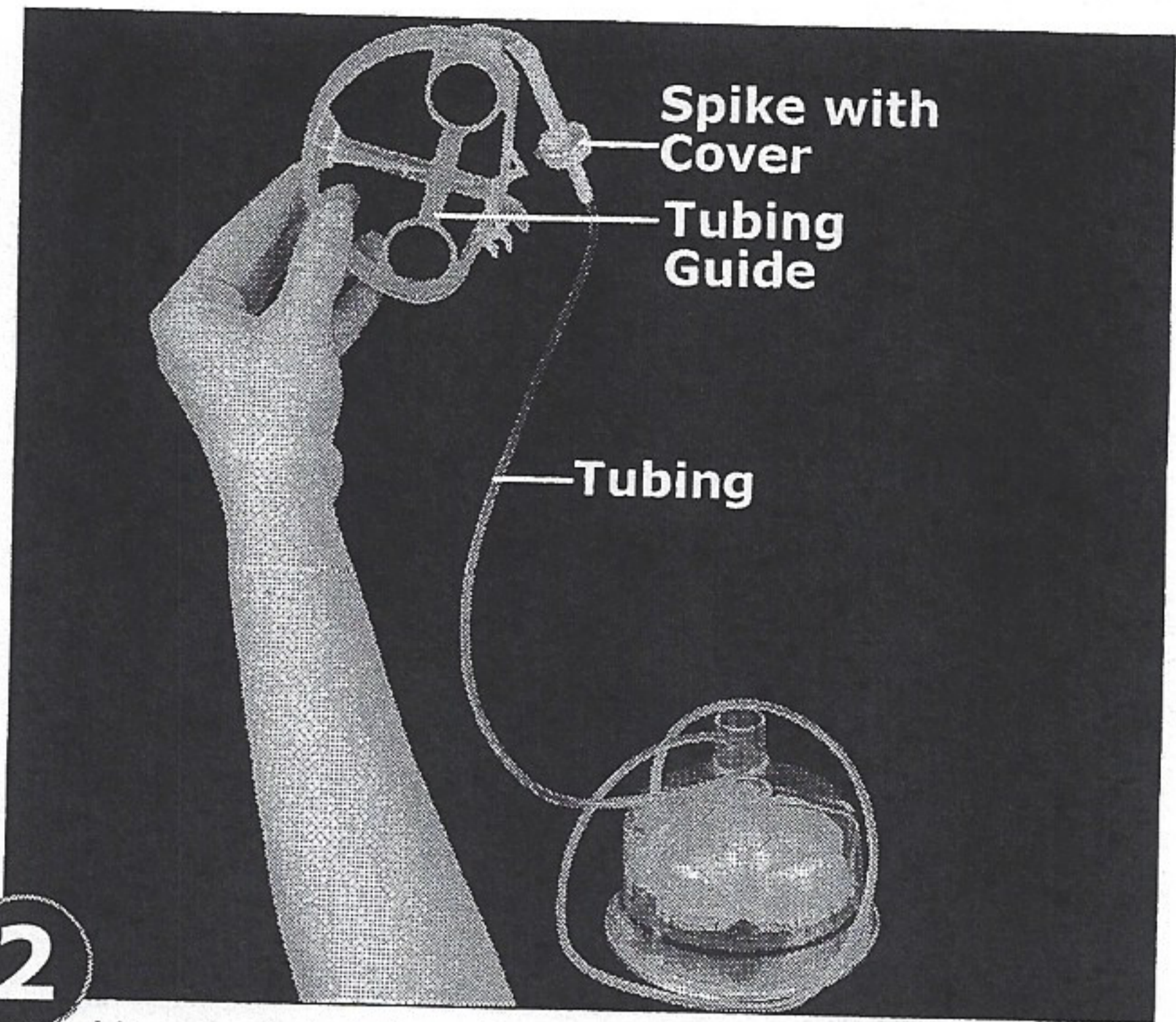

2

Unwind tubing from tubing guide.


Keep spike cover on spike until ready to use.

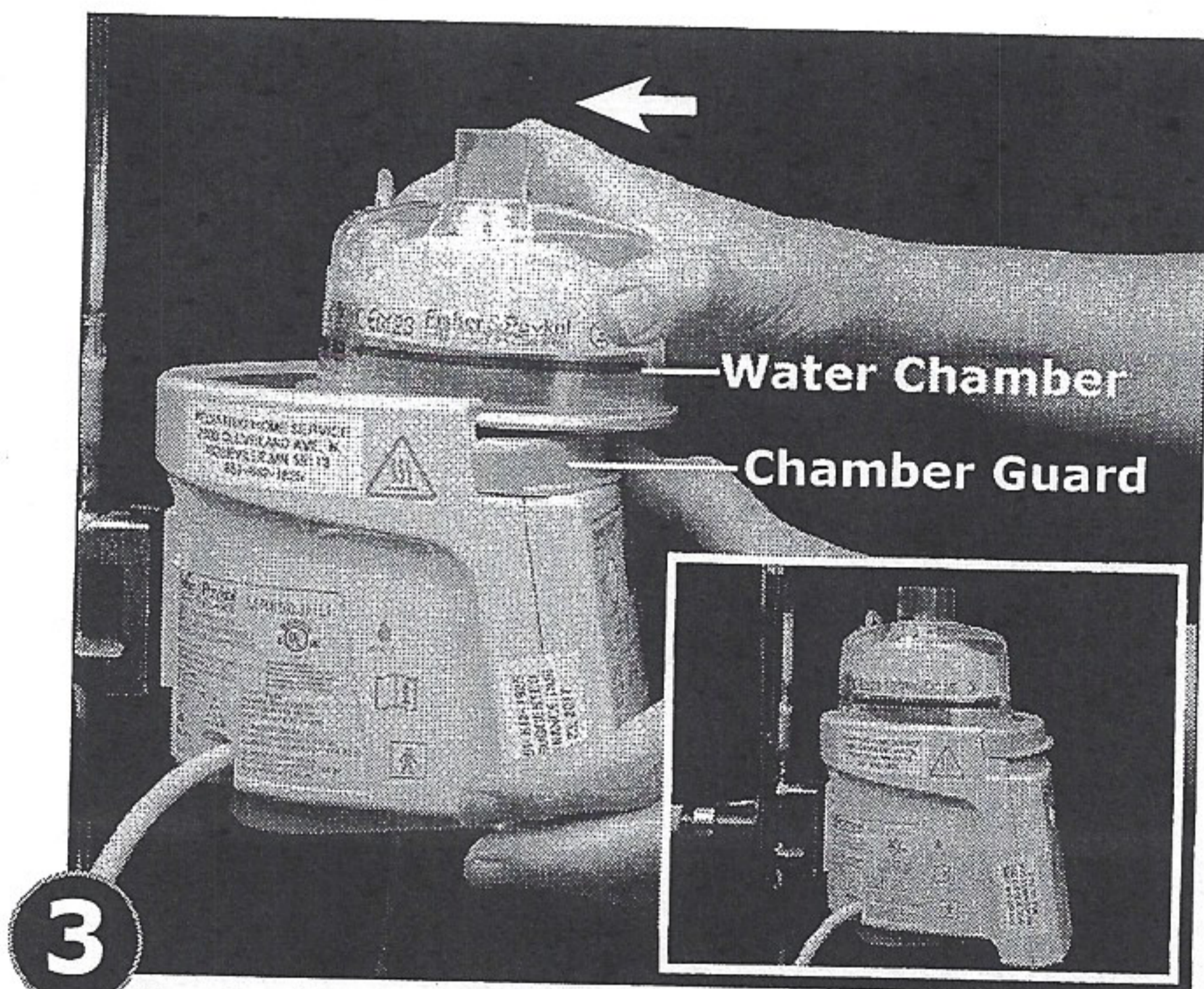

3

Push down on humidifier chamber guard.

Slide water chamber toward back of humidifier until chamber guard pops back up.

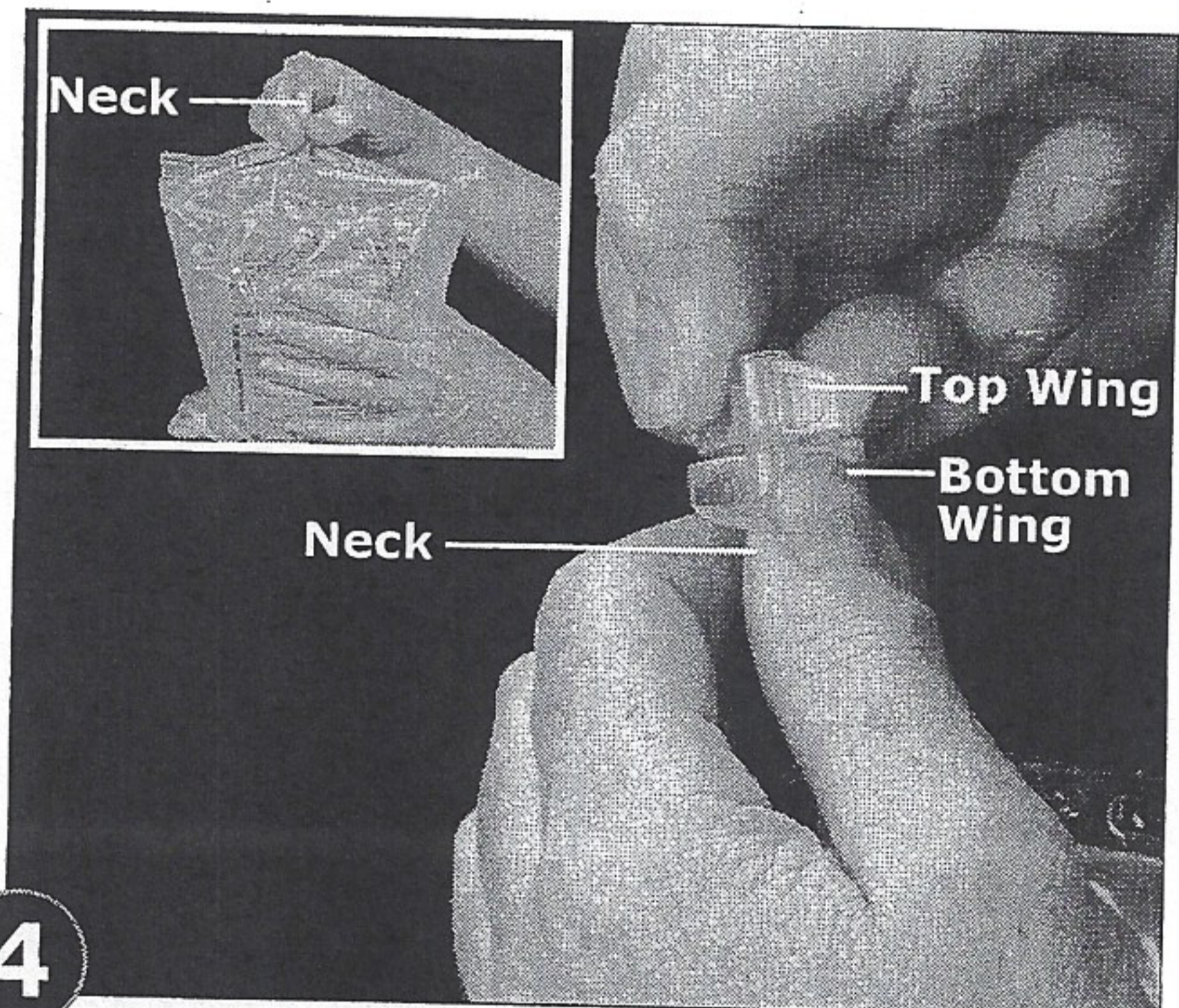

4


Hold sterile water drip bag so neck of bag points toward ceiling.

Hold neck and bottom wing on sterile water drip bag with one hand. With other hand, turn and remove top wing.

Continued on next page →


Setting Up the Humidifier and Water Chamber

[continued]

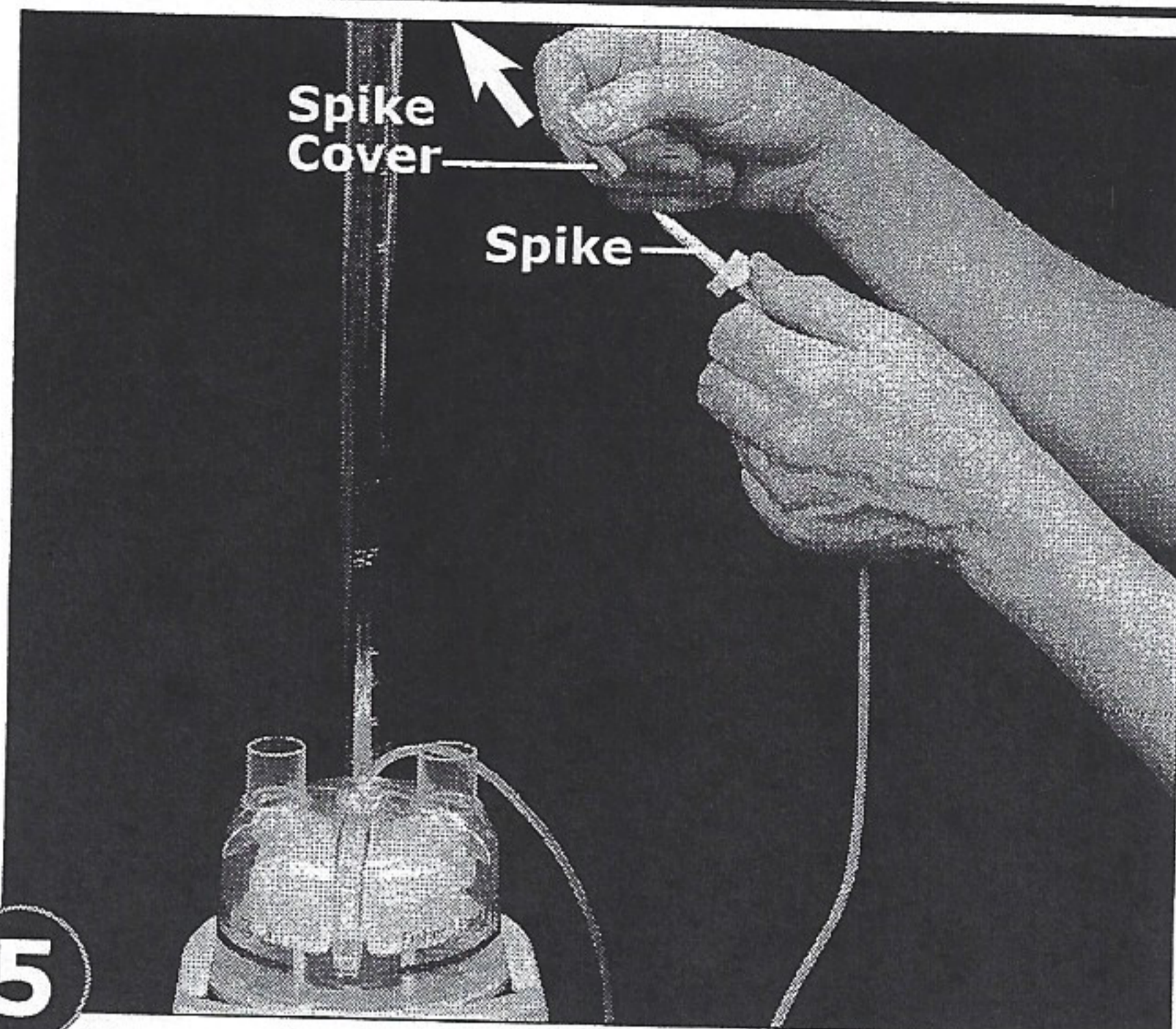


Remove cover from spike at end of water chamber tubing.

IMPORTANT! Do NOT allow spike to touch anything. *If spike touches anything, wipe with alcohol wipe.*

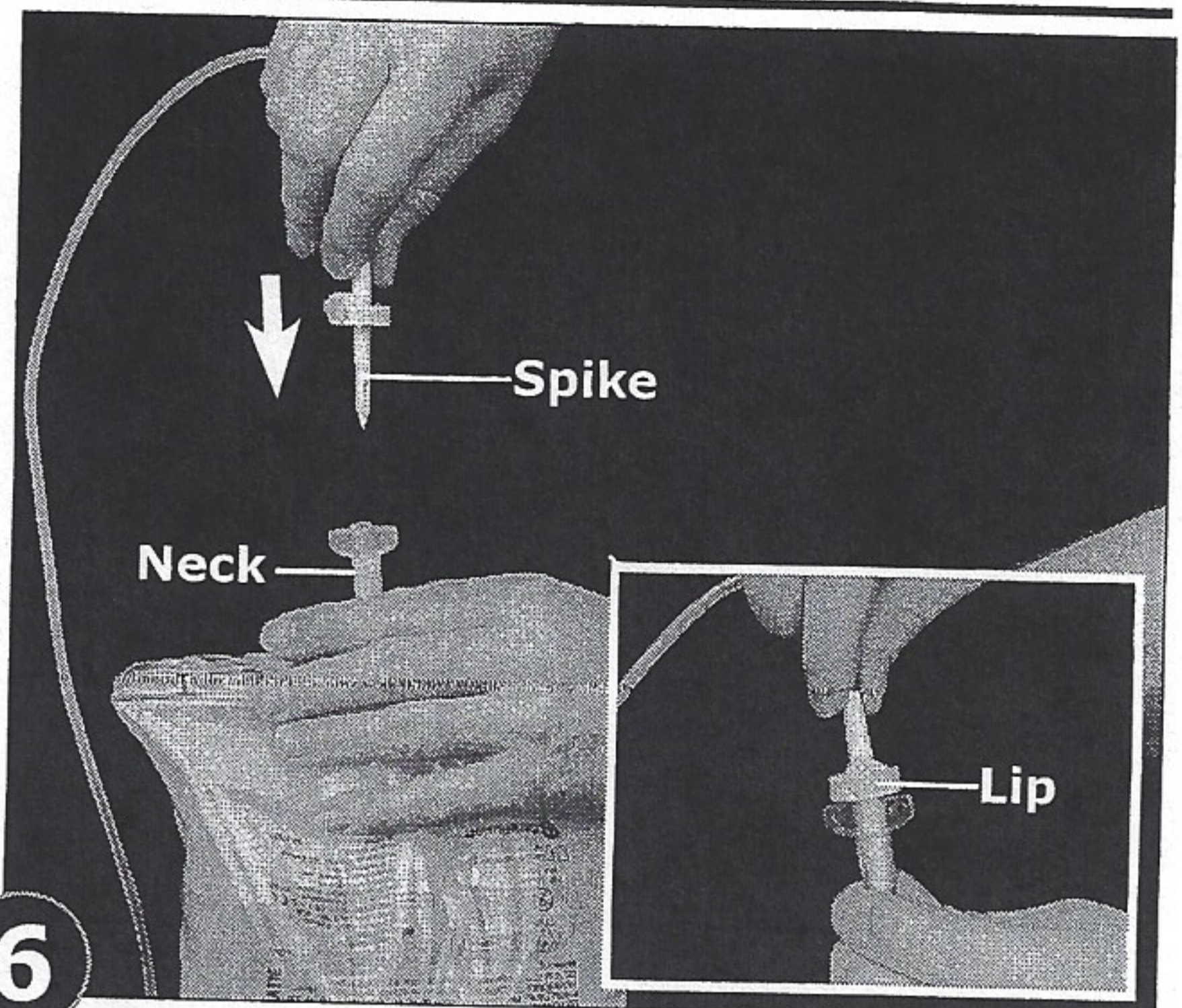

Push spike into neck of sterile water drip bag until lip on spike touches wing.

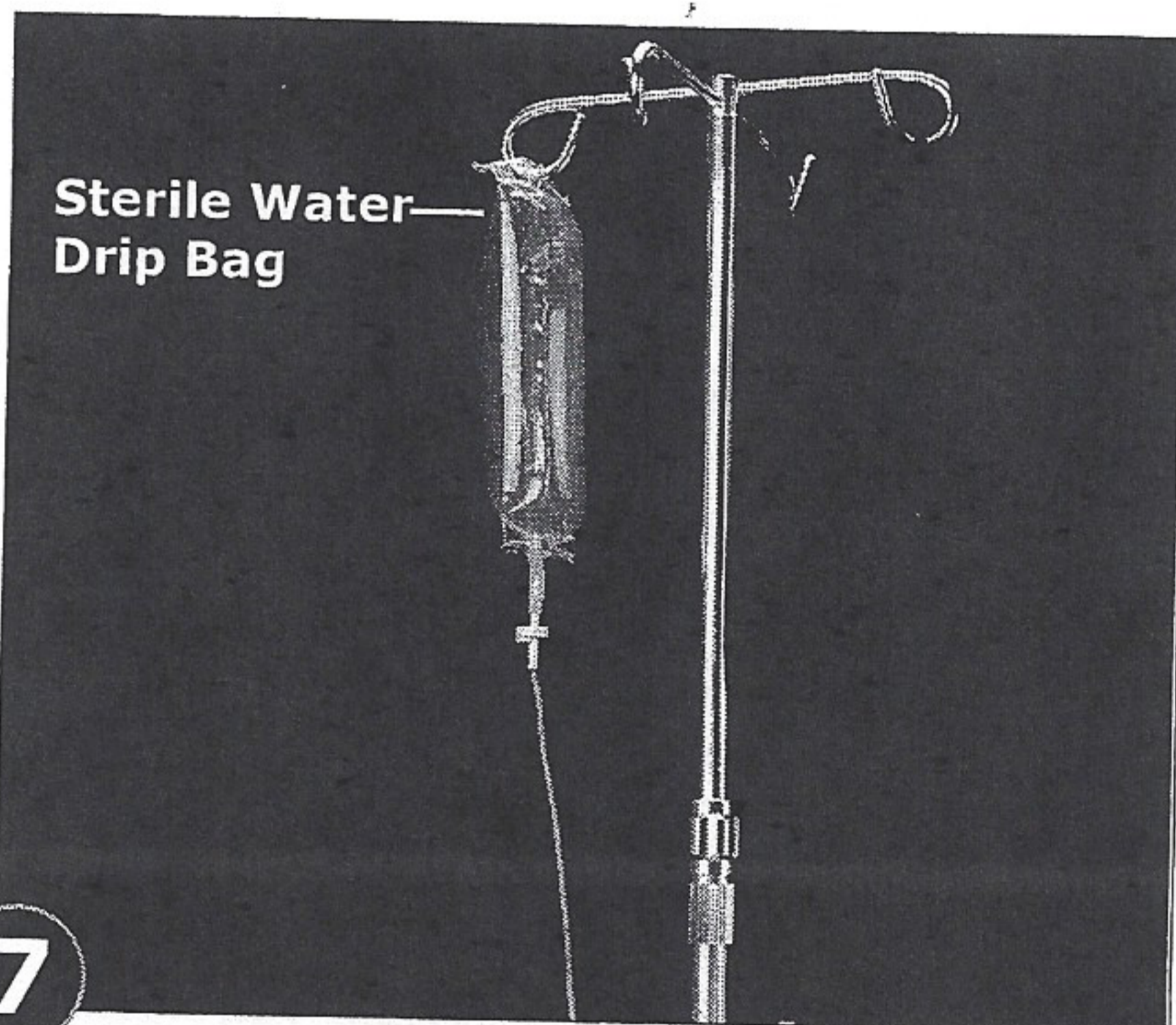

Hang sterile water drip bag from one hook on top of IV pole or stand.

Sterile water begins to fill water chamber.